

Proven & Effective - High Quality 100% Polyester Wipes with 70 | 30 IPA



Klerwipe™ 70|30 IPA

Sterile, knitted 100% polyester wipes impregnated with spore-free IPA blended with EP Water for Injection

Klerwipe 70|30 IPA has ultrasonically sealed edges and is pre-laundered to reduce the bioburden on the wipe, providing an endotoxin level of less than 1 EU per wipe.


No undesirable residues are left after using these wipes because the alcohol evaporates quickly and Klerwipe 70|30 IPA contain no chemical binders.

This makes the 100% polyester Klerwipe 70|30 IPA ideally suitable for use

on surfaces that are directly in contact with products where residues and endotoxins are unacceptable.

Klerwipe 70|30 IPA:

- ▲ Is a low residue product, no chemical binders
- ▲ Is made of heavy high-quality, weight knitted polyester wipe
- ▲ Provides high absorbency and consistent application
- ▲ Ensures quick access with the easy-to-use pouch, accessible even with gloves
- ▲ Provides a reduced risk of inhalation
- ▲ Is available in a choice of sizes
- ▲ Has ultrasonically cut and sealed wipes
- ▲ Has a lint-free particle and fibre specification

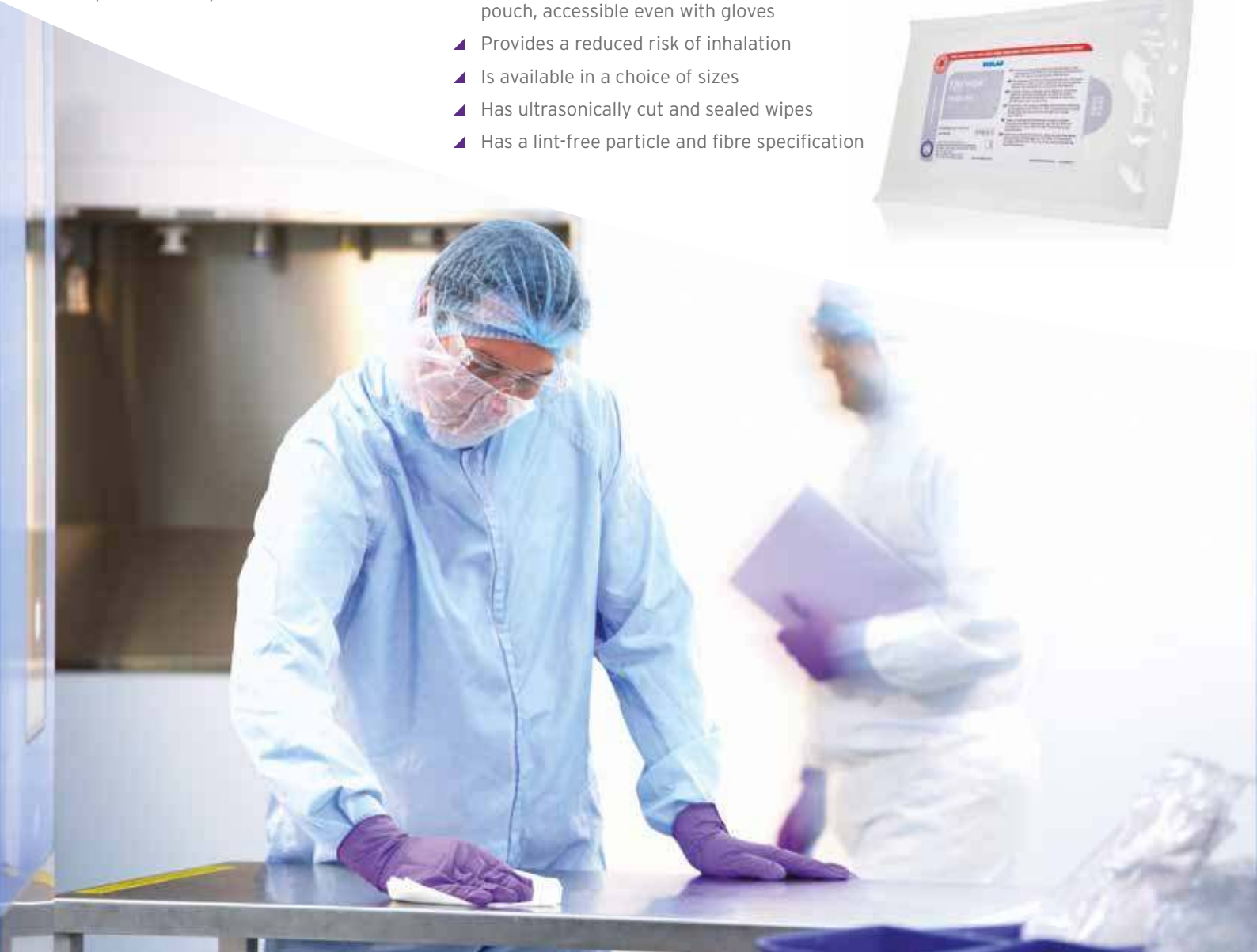



Klerwipe™ 70|30 IPA

Excellent performance low residue wipes

Why use Klerwipe 70|30 IPA?

PROVEN

Klerwipe 70|30 IPA wipes are impregnated with 70|30 v/v IPA blended with European Pharmacopoeia specification Water for Injection. They are packed and triple-bagged in a Grade C (ISO Class 7) cleanroom and irradiated using a validated process at no less than 25kGy. Each wipe has a guaranteed endotoxin level of less than 1 EU so perfect for use in critical product contact areas.

EASY & ECONOMICAL

Each pouch of Klerwipe 70|30 IPA contains just 10 low particulate wipes - ideal for sessional use without compromising sterility. The packs feature Klertear™ for easy opening, even while wearing gloves. The pouch itself does not generate any particles. Each unopened pouch has a two year shelf life from the date of manufacture.

DIRECT

Ideal for use directly on surfaces, equipment and production components - even in hard to reach areas. Klerwipe 70|30 IPA can be deployed for rapid, thorough clean down, remaining wet throughout use.

VALIDATION

Manufactured by a ISO 9001:2008, ISO 14001:2015 accredited company.

Certificates of Analysis and Irradiation are available for every batch.

Product has a two year shelf life.



Klerwipe 70|30 IPA
200 x 200mm wipe



Klerwipe 70|30 IPA
300 x 300mm wipe

MICROBIOLOGICAL AND VIROLOGICAL TEST RESULTS

EFFICACY UNDILUTED	CONTACT TIME
EN 1276 Bactericidal	5 min
EN 1650 Fungicidal	15 min
EN 13697 Bactericidal and Fungicidal	5 min

WORLDWIDE HEADQUARTERS

1 Ecolab Place
St. Paul, MN 55102
USA

www.ecolab.com

EUROPE HEADQUARTERS

Richtistr. 7
8304 Wallisellen
Switzerland

www.ecolablifesciences.com

PROVEN.
CONVENIENT.
DIRECT.

APPLICATION AREA

For use on floors, walls and hard surfaces.

INSTRUCTIONS FOR USE

For optimum results ensure complete wetting of the surfaces (max. 40 ml/m²). Re-seal when not in use. Dispose of packaging only if emptied and closed.

INGREDIENTS

100ml of liquid contains: 70 ml isopropanol in Water for Injection (WFI).

PACKAGING INFORMATION

UNIT SIZE	UNIT OF SALE
200 x 200mm Wipe	20 packs x 10 wipes
300 x 300mm Wipe	15 packs x 10 wipes

FOR FULL PRODUCT ORDER CODES, PLEASE SEE OUR PRODUCT CATALOGUE OR VISIT OUR WEBSITE

USE BIOCIDES SAFELY. ALWAYS READ THE LABEL AND PRODUCT INFORMATION BEFORE USE.

