

Kimtech™ M3 Certified Sterile Pleat-Style, Face Mask

Product Information

Code	Description	Color	Style	Case Count
62467	M3 Certified Sterile Face Mask, 7" (18cm) (Soft Ties)	White	Pleated	20/bag; 200/case
62470	M3 Certified Sterile Face Mask Knitted, 7" (18cm) (Earloops)	White	Pleated	20/bag; 200/case
62494	M3 Certified Sterile Face Mask, 9" (23cm) (Soft Ties)	White	Pleated	20/bag; 200/case

Physical Properties	Targets	Test Method
Particle Filtration Efficiency @ 0.1 micron	90%* (Avg.)	ASTM F2299
Differential Pressure in mm H ₂ O @ 8 LPM flow rate	2.4* (Avg.)	MIL-M-36954C
Bacterial Filtration Efficiency @ 3.0 micron	92%* (Avg.)	ASTM F2101
Particle Shedding (Helmke Drum)	Category 1	IEST RP-CC003.4


*Pre-sterile



Mask Construction	
Binding, Ties	White, Spunlaced Polyester
Nosepiece	Fully enclosed, soft, malleable aluminum
Outer Facing	White, Polypropylene, Apertured Film
Inner Facing	White Spunbond Fabric
Filter Media	White, Meltblown Polyester
Construction	Ultra Sonic Bonding

Additional Information:

- Manufactured in ISO 9001:2015 registered facilities.
- All Face Masks are individually bagged, then doubled bagged.at 20 aggregate.
- Recommended for cleanroom.
- Gamma irradiated with Sterility Assurance Level of 10⁻⁶.
- Certificates of Conformance available by request, send email to kcpcerts@kcc.com

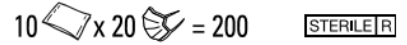

Workmanship Requirements:

Face Masks shall be made in accordance with reasonable industry practice with respect to holes, tears, wrinkles, cleanliness, odors and foreign materials or dirt.

For Additional information: Contact your Kimberly-Clark representative, our Customer Service Team (888)-346-GOKC (4652) or kimtech@kcc.com

Rev. 8/2018

Kimberly-Clark warrants that its products (1) comply with K-C's standard specifications as of the delivery date to K-C's authorized distributors/direct purchasers (2) comply with all K-C labeling representations; and (3) are manufactured in compliance with all applicable federal, state, and local laws in effect at the time and place of manufacture of the products. **THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.** K-C is not liable for any kind of special, incidental, or consequential damages. K-C's liability for breach of contract, tort or other cause of action shall not exceed the product purchase price. Purchasers and users are deemed to have accepted the above warranty and limitation of liability, and cannot change the terms by verbal agreement or by any writing not signed by K-C. To the extent required by applicable law, K-C does not limit its liability for death/injury resulting from K-C's negligence.



Public